


DuPont™ Tyvek® IsoClean® , *Model IC 105 S WH CS*

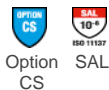

Technical Data Sheet

DuPont Tyvek® IsoClean® hooded coverall with attached slip-retardant shoe covers, model IC 105 S WH CS. Clean-processed and gamma-sterilized. Serged internal seams. Elastication at wrists, ankles and face. Thumb loops. Zipper closure. Aseptically folded. White.

Sterilisation

- This garment is not CE certified as chemical protective clothing

Packaging(Quantity/Box)

25 per box, individually packed. 2 polyethylene liners. Cardboard box.

Product Size	Article Number	Body Height(cm)	Chest Girth(cm)	Chest Girth(in)	Body Height(ft/in)
MD	D14247979	165-178	91-97	34 3/4 38 1/4	5'3"-5'7"
LG	D14247987	173-183	97-102	38 3/4 42 1/4	5'5"-5'9"
XL	D14247990	180-191	102-109	42 3/4-46 1/4	5'8"-6'2"
2X	D14248003	185-196	109-117	46 3/4-50 1/4	6'0"-6'4"
3X	D14248016	188-198	117-132	54 3/4-58 1/4	6'2"-6'4"

Reference Number: IC0105SWHCS

PHYSICAL PROPERTIES

Property	Test Method	Result	EN Class
Abrasion Resistance ⁷	EN 530 Method 2	>10 cycles	1 of 6 ¹
Basis Weight	DIN EN ISO 536	45 g/m ²	N/A
Colour	N/A	White	N/A
Exposure to high Temperature	N/A	Melting point ~135 °C	N/A
Flex Cracking Resistance ⁷	EN ISO 7854 Method B	>100000 cycles	6 of 6 ¹
Puncture Resistance	EN 863	>5 N	1 of 6 ¹
Resistance to Water Penetration	DIN EN 20811	7 kPa	N/A
Surface Resistance at RH 25%, inside ⁷	EN 1149-1	2 ¹⁰ Ohm	N/A
Tensile Strength (MD)	DIN EN ISO 13934-1	>30 N	1 of 6 ¹
Tensile Strength (XD)	DIN EN ISO 13934-1	>30 N	1 of 6 ¹
Thickness	DIN EN ISO 534	185 µm	N/A
Trapezoidal Tear Resistance (MD)	EN ISO 9073-4	>10 N	1 of 6 ¹
Trapezoidal Tear Resistance (XD)	EN ISO 9073-4	>10 N	1 of 6 ¹

¹ According to EN 14325 ² According to EN 14126 ³ According to EN 1073-2 ⁴ According to EN 14116 ¹² According to EN 11612 ⁵ Front Tyvek ® / Back
 Instructions for Use for further information, limitations and warnings > Larger than < Smaller than **N/A** Not Applicable **STD DEV** Standard Deviation ⁶ Based on test according to ASTM D-572 ⁷ See

COMFORT

Property	Test Method	Result	EN Class
Air Permeability (Gurley method)	ISO 5636-5	Yes	N/A
Air Permeability (Gurley method)	ISO 5636-5	4 s	N/A
Thermal Resistance, R _{ct}	EN 31092/ISO 11092	10*10 ⁻³ m ² *K/W	N/A
Thermal Resistance, clo value	EN 31092/ISO 11092	0.065 clo	N/A
Water Vapour Resistance, Ret	EN 31092/ISO 11092	6.8 m ² *Pa/W	N/A

2 According to EN 14126 5 Front Tyvek ® / Back > Larger than < Smaller than **N/A** Not Applicable

PENETRATION AND REPELLENCY

Property	Test Method	Result	EN Class
Repellency to Liquids, Sodium Hydroxide (10%)	EN ISO 6530	>90 %	2 of 3 ¹
Repellency to Liquids, Sulphuric Acid (30%)	EN ISO 6530	>95 %	3 of 3 ¹
Resistance to Penetration by Liquids, Sodium Hydroxide (10%)	EN ISO 6530	<5 %	2 of 3 ¹
Resistance to Penetration by Liquids, Sulphuric Acid (30%)	EN ISO 6530	<1 %	3 of 3 ¹

1 According to EN 14325 > Larger than < Smaller than

BIOLOGICAL BARRIER

Property	Test Method	Result	EN Class
Resistance to Penetration by Biologically Contaminated Aerosols	ISO/DIS 22611	Pass	1 of 3 ²
Resistance to Penetration by Blood and Body Fluids using Synthetic Blood	ISO 16603	Pass	3 of 6 ²
Resistance to Penetration by Blood-borne Pathogens using Bacteriophage Phi-X174	ISO 16604 Procedure C	No classification	No classification ²
Resistance to Penetration by Contaminated Liquids	EN ISO 22610	Pass	1 of 6 ²
Resistance to Penetration by Contaminated Solid Particles	ISO 22612	Pass	1 of 3 ²

2 According to EN 14126 > Larger than < Smaller than

CLEANILESS

Property	Test Method	Result	EN Class
Bacterial Filtration Efficiency (3 µm)	ASTM F2101	98.4 % ± 0.9 % STD DEV	N/A
Particle Shedding (Helmke Drum)	IEST-RP-CC003.4.	Category I	N/A

5 Front Tyvek ® / Back > Larger than < Smaller than **N/A** Not Applicable **STD DEV** Standard Deviation

Permeation Data for Tyvek® IsoClean®

Hazard / Chemical Name	Physical State	CAS	BT Act	BT 0.1	BT 1.0	EN	SSPR	MDPR	Cum Time ISO 480 150
Carboplatin (10 mg/ml)	Liquid	441575-94-4	>240	>240	>240	5	<0.001	0.001	
Carmustine (3.3 mg/ml, 10 % Ethanol)	Liquid	154-93-8	imm	imm	>240	5	<0.3	0.001	
Cisplatin (1 mg/ml)	Liquid	15663-27-1	>240	>240	>240	5	<0.001	0.001	
Cyclo phosphamide (20 mg/ml)	Liquid	50-18-0	>240	>240	>240	5	<0.008	0.008	
Doxorubicin HCl (2 mg/ml)	Liquid	25136-40-9	>240	>240	>240	5	<0.001	0.001	
Etoposide (Toposar®, Teva) (20 mg/ml, 33.2 % (v/v) Ethanol)	Liquid	33419-42-0	>240	>240	>240	5	<0.01	<0.01	
Fluorouracil, 5- (50 mg/ml)	Liquid	51-21-8	imm	imm	imm		na	0.001	
Gemcitabine (38 mg/ml)	Liquid	95058-81-4	imm	>60	>240	5	<0.4	0.005	
Ifosfamide (50 mg/ml)	Liquid	3778-73-2	>240	>240	>240	5	<0.009	0.009	
Oxaliplatin (5 mg/ml)	Liquid	63121-00-6	imm	imm	imm		na	0.001	
Paclitaxel (Hospira) (6 mg/ml, 49.7 % (v/v) Ethanol)	Liquid	33069-62-4	>240	>240	>240	5	<0.01	<0.01	
Thiotepa (10 mg/ml)	Liquid	52-24-4	imm	imm	imm		na	0.001	

BT Act (Actual) Breakthrough time at MDPR [mins] **BT 0.1** Normalized breakthrough time at 0.1 µg/cm²/min [mins] **BT 1.0** Normalized breakthrough time at 1.0 µg/cm²/min [mins] **EN** Classification according to EN 14325
SSPR Steady state permeation rate [µg/cm²/min] **MDPR** Minimum detectable permeation rate [µg/cm²/min] **CUM 480** Cumulative permeation mass after 480 mins [µg/cm²] **Time 150** Time to reach cumulative permeation mass of 150 µg/cm² [mins] **ISO** Classification according to ISO 16602 **CAS** Chemical abstracts service registry number **mins** Minutes **>** Larger than **<** Smaller than **imm** Immediate (< 4 min) **nm** Not tested
sat Saturated solution **N/A** Not Applicable * Based on lowest single value **na** Not attained **8** Actual breakthrough time; normalized breakthrough time is not available

Important Note

The permeation data published have been generated for DuPont by independent accredited testing laboratories according to the test method applicable at that time (EN ISO 6529 (method A and B), ASTM F739, ASTM F1383, ASTM D6978, EN369, EN 374-3)

The data is typically the average of three fabrics samples tested.

All chemicals have been tested at an assay of greater than 95 (w/w) % unless otherwise stated.

The tests were performed between 20 °C and 27°C and at environmental pressure unless otherwise stated.

A different temperature may have significant influence on the breakthrough time.

Permeation typically increases with temperature.

Cumulative permeation data have been measured or have been calculated based on minimum detectable permeation rate.

Cytostatic drugs testing has been performed at a test temperature of 27°C according to ASTM D6978 or ISO 6529 with the additional requirement of reporting a normalized breakthrough time at 0.01 µg/cm²/min.

Chemical warfare agents (Lewisite, Sarin, Soman, Mustard, Tabun and VX Nerve Agent) have been tested according to MIL-STD-282 at 22°C or according to FINABEL 0.7 at 37°C.

Permeation data for Tyvek® is applicable to white Tyvek® 500 and Tyvek® 600 only and is not applicable for other Tyvek® styles or colours.

Permeation data are usually measured for single chemicals. The permeation characteristics of mixtures can often deviate considerably from the behaviour of the individual chemicals.

The permeation data for gloves published have been generated according to ASTM F739 and to ASTM F1383.

The degradation data for gloves published have been generated based on a gravimetric method.

This degradation testing exposes one side of the glove material to the test chemical for four hours. The percent weight change after exposure is measured at four time intervals: 5, 30, 60 and 240 minutes.

Degradation Ratings:

- E: EXCELLENT (0-10% Weight Change)
This garment does not protect against ionizing radiation.
- D: GOOD (11-20% Weight Change)
This garment is not flame resistant and should not be used around heat, open flame, sparks or in potentially flammable environments.
- C: FAIR (21-30% Weight Change)
- B: POOR (31-40% Weight Change)
- A: POOR (41-50% Weight Change)
- F: FAIR (21-30% Weight Change)

The information provided herein corresponds to our knowledge on the subject at the date of its publication. This information may be subject to revision as new knowledge and experience becomes available. The data provided fall within the normal range of product properties and relate only to the specific material designated; these data may not be valid for such material used in combination with any other materials or additives or in any process, unless expressly indicated otherwise. The data provided should not be used to establish specification limits or used alone as the basis of design; they are not intended to substitute for any testing you may need to conduct to determine for yourself the suitability of a specific material for your particular purposes. Since DuPont cannot anticipate all variations in actual end-use conditions DuPont makes no warranties and assumes no liability in connection with any use of this information. Nothing in this publication is to be considered as a license to operate under or a recommendation to infringe any patent rights.

Please use the permeation data provided as a part of the risk assessment to assist with the selection of a protective fabric, garment, glove or accessory suitable for your application.

Breakthrough time is not the same as safe wear time. Breakthrough times are indicative of the barrier performance, but results can vary between the test methods and laboratories. Breakthrough time alone is insufficient to determine how long a garment may be worn once the garment has been contaminated. Safe user wear time may be longer or shorter than the breakthrough time depending on the permeation behaviour of the substance, the toxicity of the substance, working conditions and the exposure conditions (e.g. temperature, pressure, concentration, physical state).

Latest Update Permeation Data: 15/03/2019

The information provided herein corresponds to our knowledge on the subject at the date of its publication. This information may be subject to revision as new knowledge and experience becomes available. The data provided fall within the normal range of product properties and relate only to the specific material designated; these data may not be valid for such material used in combination with any other materials or additives or in any process, unless expressly indicated otherwise. The data provided should not be used to establish specification limits or used alone as the basis of design; they are not intended to substitute for any testing you may need to conduct to determine for yourself the suitability of a specific material for your particular purposes. Since DuPont cannot anticipate all variations in actual end-use conditions DuPont makes no warranties and assumes no liability in connection with any use of this information. Nothing in this publication is to be considered as a license to operate under or a recommendation to infringe any patent rights.

For further product information, literature and as well as assistance in locating a local supplier, please visit:

www.safespec.dupont.co.uk

DuPont Personal Protection

DuPont de Nemours (Luxembourg) S.à.r.l.

L-2984 Luxembourg

Tel.: +800 3666 6666 (international toll-free)

Fax: +352 3666 5071

E-mail: personal.protection@lux.dupont.com

The footnotes can be found on the SafeSPEC® website.

Copyright © DuPont. All rights reserved. The DuPont Oval Logo, DuPont®, The miracles of science® and all products denoted with ® or ® are registered trademarks or trademarks of E. I. du Pont de Nemours and Company or its affiliates.